Certain® Screw Removal Tool Kit (ISRT10)

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Indications:

In the unlikely event of screw fracture, the Certain® Screw Removal Tool Kit (ISRT10) facilitates the removal of a broken screw from the internal threads of the Certain Implant.

Contraindications:

The Certain Screw Removal Tool Kit is not to be used to remove screws from an External Hex Implant.

Warnings:

The Certain Screw Removal Tool Kit and its components are intended to be used only with BIOMET 3i Certain Implant screws.

Precautions:

Slight to moderate pressure is to be applied when using the components.

Reuse of BIOMET 3i products that are labeled for single-use may result in product contamination, patient infection and/or failure of the device to perform as intended.

Cleaning and Sterilization:

Since the components of the Certain Screw Removal Tool Kit are reusable, these must be cleaned and sterilized prior to use. Refer to P-IFSCSS for complete information.

Side Effects:

When used according to these instructions, there are no known side effects.

Caution:

Federal law restricts this device to sale by or on the order of a licensed dentist.

Components of the ISRT10 are:

ISRT01
Certain Hand
Tap

ISRT05
Certain Guide
Handle


ISRT06
Screw Removal
Tool


ISRT07
Manual Reversing
Drill


IWSU30
Certain Waxing
Screw/Guide Pin

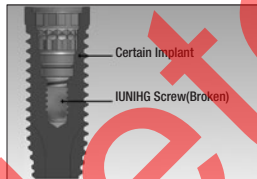**Instructions:**

Step 1: Remove abutment and coronal end of broken screw from implant.

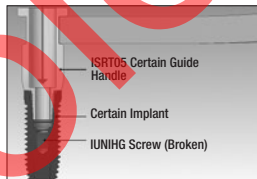
Step 2: Align and insert the ISRT05 Certain Guide Handle into the broken screw site. Seat the handle assembly fully on the platform of the implant with the hex aligned. Note: The ISRT05 Guide Handle fits both standard and 3.25mm(D) Certain Implants.

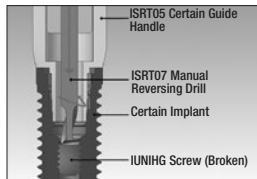
Step 3: Insert the ISRT07 Manual Reversing Drill into the ISRT05 Guide Handle. Firmly seat the drill point(s) onto the surface of the broken screw inside the implant. Slowly (with slight to moderate pressure on the screw) turn the ISRT07 Manual Reversing Drill one to two revolutions in the reverse direction (counterclockwise). Do not continue to loosen the screw from the implant.


Step 4: Remove the ISRT07 Manual

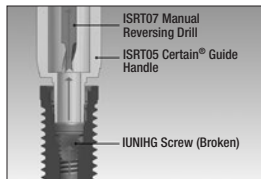Reversing Drill and then remove the ISRT05 Guide Handle from the implant.

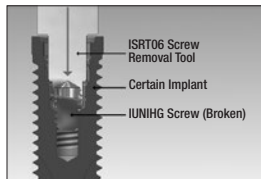
Step 5: Insert the ISRT06 Screw Removal Tool into the implant and seat the tool onto the top portion of the broken screw. Press onto the screw with slight to moderate pressure to engage/capture the screw. Rotate the tool in the reverse direction (counter-clockwise) several rotations. The screw should release from the implant internal threads.

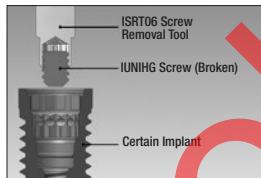
Step 6: Remove the ISRT06 Screw Removal Tool from the implant. The fractured screw should remain captured in the ISRT06. The screw can be removed by pulling it out of the tool tip. If the tool fails to catch the screw, repeat steps 2 – 5 until the screw is extracted from the implant.

Step 7: After the broken screw has been removed, insert the IWSU30 Certain

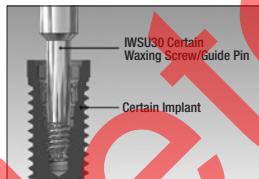Waxing Screw/Guide Pin to evaluate the integrity of the internal threads of the implant. If the IWSU30 does not rotate easily into the implant, proceed to Step 8. If the IWSU30 rotates easily into the implant, you may proceed in placement of the new abutment.

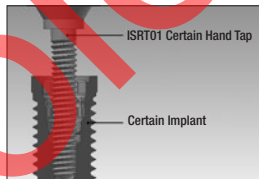Step 8: If the IWSU30 did not rotate easily into the internal threads of the implant, insert the ISRT01 Certain Hand Tap into the implant. Rotate the ISRT01 in a clockwise direction to rethread the implant. Before replacing the abutment, checking the internal threads with the IWSU30 Certain Waxing Screw/Guide Pin is recommended. Once the integrity of the implant threads has been verified, the abutment may be replaced.

LOT

Batch code

Caution, consult accompanying documents

Consult Instructions for Use
www.ifu.biomet3i.com

Use by

Legal Manufacturer

Date of Manufacture

REF

Catalogue number

By prescription only

EC REP

Authorized European Representative

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Certain is a registered trademark of BIOMET 3i LLC. ©2016 BIOMET 3i LLC.
Certain è un marchio depositato di BIOMET 3i LLC. ©2016 BIOMET 3i LLC.
Certain é uma marca comercial registrada da BIOMET 3i LLC. ©2016 BIOMET 3i LLC.
Certain est une marque déposée de BIOMET 3i LLC. ©2016 BIOMET 3i LLC.
Certain ist eine eingetragene marke von BIOMET 3i LLC. ©2016 BIOMET 3i LLC.
Certain es una marca comercial registrada de BIOMET 3i LLC. ©2016 BIOMET 3i LLC.